detax

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Precision impression silicone for dental prosthetics

INDICATION

silaplast Futur, hydro C Putty, exaplast putty – tray material for correction impressions (two-step), double-mix technique (single-step) and anatomic impressions.

silasoft® Normal, silasoft® Special, hydro C Wash, exasoft wash – correction material, sprayable, for correction impressions (two-step), double-mix technique and anatomic impression (single-step), relining impression only for silasoft® Special and silasoft® Normal.

silone® – for taking impressions of dentate jaws, anatomic impressions, onlay and inlay impressions.

sta-seal – for taking impressions of edentulous jaws, mucostatic functional and relining impressions.

sta-seal f – for functional impression-taking with extended curing phase, functional margin design.

PATIENT TARGET GROUP

Persons being treated in the context of a dental procedure.

Do not use for the treatment of pregnant/breastfeeding women and children (catalyst component contains alkyl stannosiloxanes).

INTENDED USERS

Dentist

PRODUCT DESCRIPTION

The CDS silicones are condensation-curing, elastomeric impression materials.

The CDS silicones have differing consistencies in accordance with DIN EN ISO 4823 type 0 – 3 (see table on page 16).

METERING AND MIXING (see table on page 16)

silaplast Futur – with catp universal FUTUR, paste / with silaplast Futur catf, liquid

silasoft® Normal, silasoft® Special – with cat p universal Futur, paste / with silasoft® catf, liquid

silone®, sta-seal – with catp universal FUTUR, paste / with sta-seal catf, liquid

sta-seal f – with sta-seal catf, liquid

hydro C Putty, hydro C Wash – with hydro C Activator, paste

exaplast putty, exasoft wash – with exagel cat

silaplast Futur, hydro C Putty, exaplast putty, sta-seal f – manual metering and kneading

silasoft® Normal, silasoft® Special, silone®, sta-seal, hydro C Wash, exasoft wash – manual metering and mixing with a spatula

PROCEDURE

TRAY, ADHESIVE

Select a suitable impression tray in accordance with the impression technique. To optimize silicone adhesion, apply sili Adhesive Fluid or spray it on (Fig. 1).

MANUAL METERING AND MIXING

Putty – The metering is carried out with the enclosed metering spoon (Fig. 2a) or with the dosing syringe (Fig. 2b).

Paste – The metering is carried out from the tube onto a mixing block in the required strand length (Fig. 3).

Catalyst liquid – The metering is carried out with the enclosed metering spoon or by dripping, respectively.

Mixing is done manually by kneading or using a spatula on a mixing block until a homogeneous mixture or uniform color, respectively, is achieved. The mixing time is 30 seconds (Fig. 4). **silasoft® N / silasoft® S / hydro CWash / exasoft wash** can be taken up and applied using a disposable or an impression syringe (Fig. 5a/5b). For **sta-seal f**, use of an easy Application Syringe is recommended for individual margin design (Fig. 5c). After metering, immediately close the container again.

DISINFECTION

After taking the impression out of the mouth, rinse it under running lukewarm water. Subsequent disinfection can be carried out with suitable commercially available disinfecting solutions in accordance with the manufacturer's instructions.

PRODUCING THE LABORATORY MODEL

With dental plasters type III – V, this can be done already about 30 minutes, up to no later than 72 h, after removal from the mouth.

TRAY CLEANING

Remove cured material mechanically. Remove sili Adhesive Fluid with sili Adhesive Solvent. Use only in well-ventilated areas. Then clean and disinfect the tray as usual.

SAFETY INFORMATION

- ▶ Only for the specified intended use by trained specialists.
- ▶ Wear personal protective equipment (protective gloves, goggles) when handling the uncured material.
- ▶ Wear suitable personal protective equipment (protective gloves, goggles, face mask) when finishing the cured material.
- ▶ Avoid eye contact with catalyst, rinse eyes immediately with plenty of flowing water if necessary, consult ophthalmologist immediately.
- ▶ Cured impression materials are chemically resistant – avoid stains on clothing.
- ▶ Do not leave any residues of the impression material in the patient's mouth or sulcus, and make sure to avoid inadvertent pressing of material into areas outside the oral cavity through mouth-antrum connections or malformations of the palate.
- ▶ Refer to the relevant safety data sheet for hazard and safety information.

NOTES

- ▶ Detax shall not be held liable for any damage caused by misuse.
- ▶ The good flowability of **silasoft® N** ensures excellent detail fidelity. Because of this very light-bodied consistency, in rare cases **silasoft® N** is prone to partial segregation. Before metering the catalyst, homogenize any segregated **silasoft® N** with a spatula on the mixing block or in the bottle.
- ▶ **exasoft wash / hydro CWash**: Any silicone oil segregated by storage can be mixed normally with the other components and has no effect on processing. The curing behavior is not affected thereby.
- ▶ Overmetering of catalyst and higher temperatures accelerate, undermetering and lower temperatures delay the curing.
- ▶ Carefully close the catalyst bottle or tube immediately after use.
- ▶ Read and understand the safety data sheet!

For users and/or patients:

Any serious incidents occurring in relation to this product should be reported immediately to incident@detax.com and to the competent authority of the Member State in which the user and/or patient is established.

CONTRAINDICATION





Do not use this material in case of known allergies to any of the ingredients, or contact allergies. Silicone impression materials have been tried and tested millions of times; undesired effects are not to be expected in case of proper application. However, as a matter of principle immune reactions such as allergies or irritations cannot be ruled out. In case of doubt, we recommend performing an allergy test before application.

SIDE EFFECTS

Product may cause allergic reactions.

DISPOSAL

Disposal of the contents/container must be carried out in accordance with the local/regional/national and international regulations.

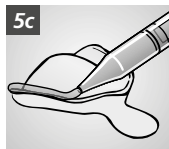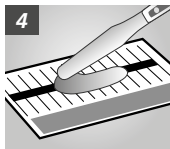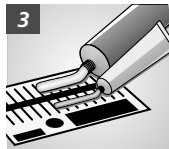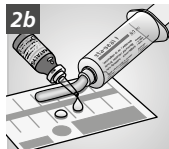
Technical Specifications	DIN EN ISO 4823 TYPE	Dosage Base mass	Dosage cat f	Dosage cat p	Product color Base/Catalyst
silaplast Futur	type 0, putty	1 ×  = 13.5 g	1 ×  = 0.5 g		white / light blue
				6 cm  = 0.3 g	white / blue
silasoft® Normal	type 3, light-bodied	10 cm  = 5 g	10 ×  = 0.2 g		pink opaque / blue
				10 cm  = 0.5 g	pink opaque / blue
silasoft® Special	type 3, light-bodied	10 cm  = 4.5 g	8 ×  = 0.16 g		pink-transp. / blue
				10 cm  = 0.5 g	pink-transp. / blue
silone®	type 2, medium-bodied	1 cm  = 1 g	1 ×  = 0.02 g		russet / blue
				1 cm  = 0.05 g	russet / blue
sta-seal	type 2, medium-bodied	1 cm  = 1.1 g	1 ×  = 0.02 g		lime green / blue
				1 cm  = 0.05 g	lime green / blue
sta-seal f	type 1, heavy-bodied	1 cm  = 1.7 g	1 ×  = 0.02 g		green / blue
hydro C Putty	type 0, putty	1 ×  = 13.5 g		8 cm  = 0.4 g	pastel green / blue
hydro C Wash	type 3, light-bodied	10 cm  = 5 g		10 cm  = 0.5 g	pink opaque / blue
exaplast putty	type 0, putty	1 ×  = 13.5 g		6 cm  = 0.3 g	yellow / blue
exasoft wash	type 3, light-bodied	10 cm  = 5 g		10 cm  = 0.5 g	pink / blue

Open time* min. ≤	Time in mouth min. approx.	Curing time* min. approx.
1:15	2:30	3:45
1:15	3:30	4:45
1:15	2:30	3:45
1:15	2:30	3:45
1:30	2:30	4:00
1:30	2:30	4:00
1:15	2:30	3:45
1:15	2:30	3:45
1:15	2:30	3:45
1:15	2:30	3:45
2:00	4:00	6:00
1:15	3:30	4:45
1:15	2:45	4:00
1:15	3:30	4:45
1:15	2:30	3:45


Processing: At $23^{\circ}\text{C} \pm 2^{\circ}\text{C}$,
 $50 \pm 5\%$ relative humidity

Storage:

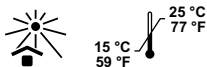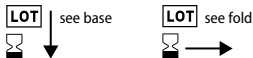* from start of mixing at $23^{\circ}\text{C} \pm 2^{\circ}\text{C}$,
 $50 \pm 5\%$ rel. humidity. Higher temperatures shorten, lower temperatures extend the indicated times.

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detax GmbH